


User's Guide

Thank you very much for purchasing MT-350R, high speed airturbine handpiece. This handpiece is designed for easy and simple use. You can utilize it to its fullest potential and enjoy for a long time without trouble if you use it properly.

MT-350R HIGH SPEED HANDPIECE

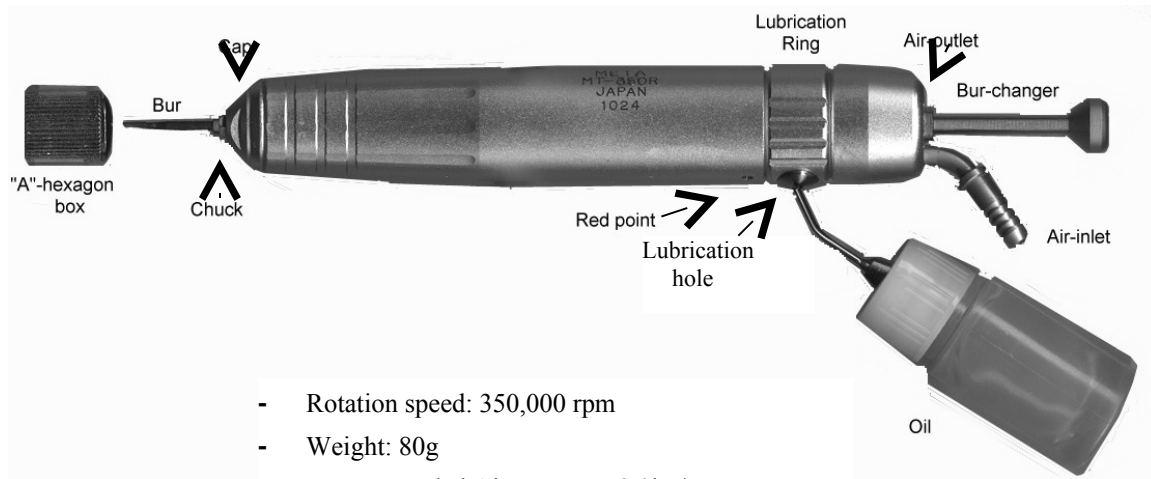

- Rotation speed: 350,000 rpm
- Weight: 80g
- Recommended Air pressure: 3.0kg/ square cm
- Burs: Friction grip(FG) 1.60 mm diameter
- Lubrication: High quality spindle oil

Inserting and removing the bur

- Inserting: Insert the bur into the chuck, and place the head of the bur in the small hole in the round side of "A" hexagon box and push the bur into the chuck side.
- Remove: Push the bur-changer
- The failure of insertion the bur to the end could cause its jumping-out


Frequency of Lubrication

- A few drops into the Lubrication hole every 2 hour use
- Use high quality spindle oil


Distributed by

How to change cartridge and chuck


- ① Loosen and take off the cap by wrench.

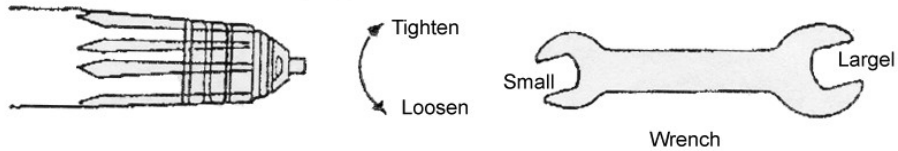


- ② Pull off the cartridge.

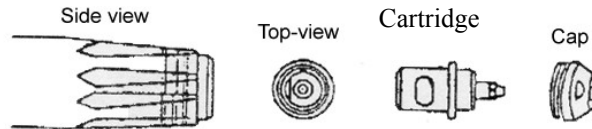


- ③ Join A-hexagon box, cartridge, B-driver, together. On squeezing all together, take off the chuck-holder by rotating B-driver to clock-wise.

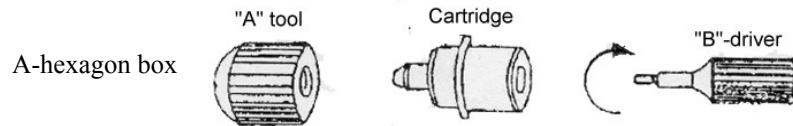

- ④ Push off the chuck with C-chuck releaser.

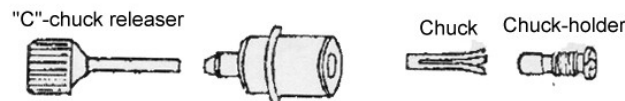

- ⑤ Install the new chuck into cartridge with reversed-order operation of the above.

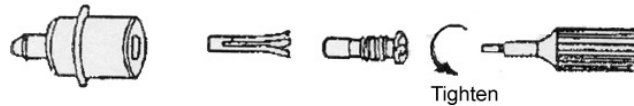

- ⑥ Push to the end the projected round-part on cartridge along the ditch of inside-body. **Very Important**

